

WOMENS HEALTH



SILICONE PESSARIES



Ring

The Ring pessary, available with or without a supporting membrane, is primarily used for mild first-degree uterine prolapse and cystocele.



Donut

The Donut pessary is primarily used for second or third-degree uterine prolapse, cystocele, and rectocele. This pessary, often used by post-menopausal women, is easily compressed for insertion and removal.



Cube

The Cube pessary, available with or without drains, is primarily used for cystocele, rectocele, second or third-degree uterine prolapse, and complete procidentia.



Gellhorn


The Gellhorn pessary with drains is primarily used for cystocele, rectocele, and second or third-degree uterine prolapse. This pessary's stem with knob is easily folded to aid with insertion and removal.

MVA Plus® Aspirator



Ipas MVA Plus® Aspirator

For uterine evacuation Product Code: PLUS

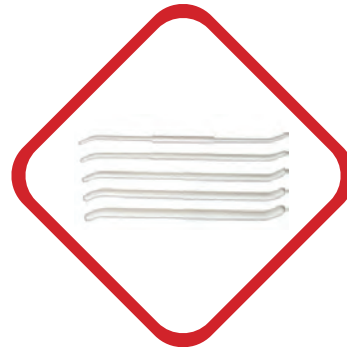

Denniston Dilators

For cervical dilation associated with uterine evacuation Product Code: DD5-S



Karman Cannulae

Flexible Karman Cannulae

8 cannulae sizes 4-10mm and 12mm. Sizes 4-6mm directly compatible with the Ipas Single Valve Aspirator. All Flexible Karman cannulae require adapters for use with Ipas MVA Plus® or Ipas Double Valve Aspirator. Disposable and packed sterile. Softer, more flexible cannula may be better option for retroverted or retroflexed uterine position.



Cannulae Adaptor Sets

For use in conjunction with applicable cannulae and MVA Aspirators
Product Codes: ADSET (sizes 6-10 MM)
ADP6MM-ADP 10MM (set of 5)



Ipas EasyGrip® Cannulae

Semi-rigid
Ready to use and sterile packed
Permanently integrated winged base,
color-coded by size
Printed insertion guides
Two opposing apertures (4-8 mm)

Large, single scoop aperture (9, 10, 12 mm)
Approximately 24 cm (9.5 in) in length
Available in eight diameters: 4, 5, 6, 7, 8, 9, 10, and 12 mm
Latex-free plastic
FDA-listed and compliant with ISO 13485 and CE-marked

IPAS EASYGRIP AND FLEXIBLE KARMEN CANNULAE



Ipas MVA Plus® Aspirator

For uterine evacuation

EFFECTIVE · SAFE · AFFORDABLE · DISCREET



MANUAL VACUUM ASPIRATION

Clinical Indications

- Treatment of incomplete abortion for uterine sizes up to 12 weeks since the woman's last menstrual period

Discreet and Patient-Friendly

- Quiet procedure meets patients' needs for privacy
- Enables faster recovery times compared to operating room procedures

Safe and Effective

- 98% effective with low complication rates*
- Safe alternative to sharp curettage

Affordable and Versatile

- Reduces procedural costs by performing procedures in an outpatient setting
- Low initial investment and per-use cost

* In less than 2% of cases one or more of the following can occur during or after the procedure: uterine or cervical injury, pelvic infection, vaginal reaction, incomplete evacuation, or acute hematometra

IPAS MVA PLUS® ASPIRATOR

- Reusable
- Steam autoclavable at 121° C (250° F)
- Processable by cold methods such as Cidex® and Sporox® II
- Easy disassembly and reassembly
- 60 cc volume
- 610-660 mm (24-26 in) Hg vacuum
- Latex-free plastic
- FDA-listed, CE-marked, and ISO 13485 compliant



WomanCare Global Trading CIC is a UK Registered Community Interest Company (No.07231920).
Email: contact@dktwomancare.org

©WomanCare Global, 2018
A-CAT-NPS-501 Rev. 00

WomanCare GLOBAL A COMPANY OF **dkd** INTERNATIONAL

CERVICAL CANCER PROGRAMME

Cryotherapy System MGC-200

Product Information


019000 MedGyn Cryotherapy System
MGC-200

Description

The MedGyn Cryotherapy System was designed by doctors and engineered by cryogenics to ensure performance, safety, and convenience.



MedGyn's Latest Innovation in Cervical Cancer Treatment **MTA-100 Thermal Ablation System**

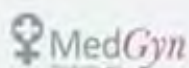

t: +1 630.627.4105 e: info@medgyn.com w: medgyn.com
Patent Pending

Description

In 2019, WHO issued new guidelines for the use of thermal ablation for cervical pre-cancerous lesions. While they continue to recommend cryotherapy, the WHO has offered thermal ablation as an alternative method of treatment for specific environments.

MedGyn is pleased to announce its latest device:

MedGyn's MTA-100 Thermal Ablation System.

The device is a portable, lightweight device for the treatment of pre-cancerous lesions optimized for 20, 40, or 60 second intervals of treatment.

ELECTROSURGICAL GENERATOR SYSTEM

Product Information

018011 MedGyn Electrosurgical Generator System

Description

The MedGyn Electrosurgical Generator System is a complete system for all LEEP/LLETZ procedures. The full system includes the Electrosurgical Generator, Smoke Evacuator and Trolley. MedGyn's Electrosurgical Generators feature both monopolar and bipolar functions to meet all your surgical demands with safety, flexibility, reliability, and convenience.



AL-101 COLPOSCOPE

Product Information

010199 MedGyn AL-101 Colposcope

Description

MedGyn's AL-101 (1-step magnification) colposcope is easy to operate, with excellent optics and the preferred option for smaller practices looking to benefit from low-cost options





DISPOSABLE COLPOSCOPY SET

Product Information

030090 Disposable Colposcopy Set

030091 Disposable Colposcopy Set with Rotating Biopsy Punch

Description

The kit contains everything needed from PVP solution to a self-lighted speculum. The speculum's integrated battery-operated LED light is cool.

STAINLESS STEEL SPECULUMS



**Weisman-Graves Vaginal Speculum
Medium**
R/Open SKU: SP201



**Graves Vaginal Speculum
Large**
SKU: SP003



**Pederson Vaginal Speculum
Large**
SKU: SP105

DISPOSABLE SPECULUMS



LLETZ SPECULUMS




**Lletz Graves Vaginal Speculum w/
Evacuation Tube
Large**




**Lletz Thumb Dress Forceps 8' SKU:
LZ101**

Contrelle® Activgard

Bladder support for women

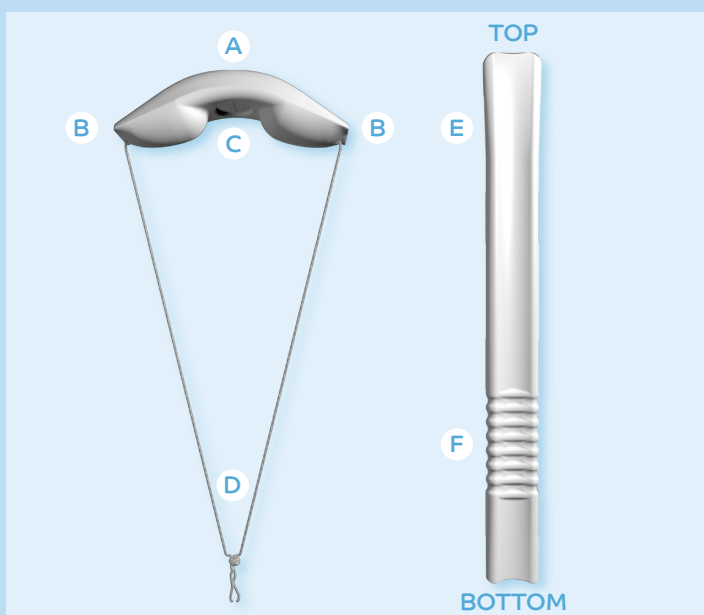

- A Contrelle® Activgard
- B Wings
- C Central depression
- D String
- E Applicator
- F Placement grooves

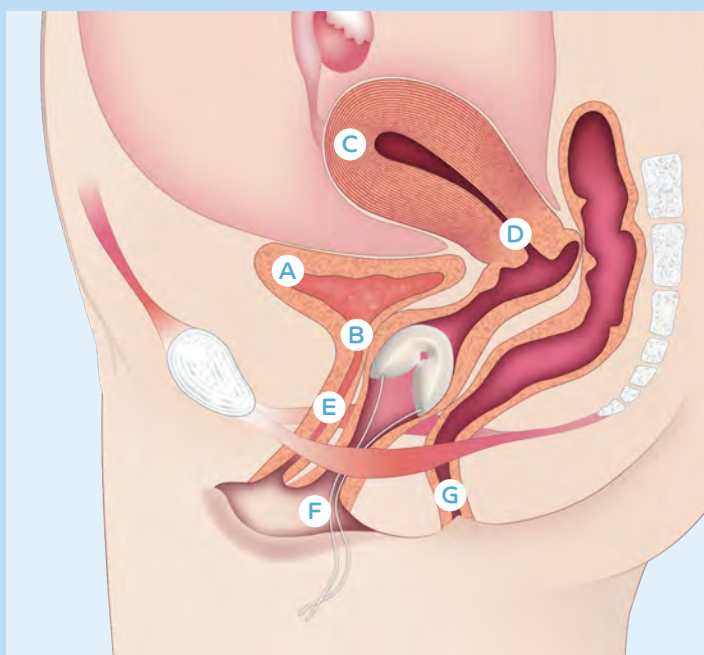

- A Bladder
- B Bladder Neck
- C Uterus
- D Cervix
- E Urethra
- F Vaginal Opening
- G Rectum



FOR ORDERS
Email: supplies@kumis.co.za
www.kumis.co.za